


EN 08 Perforated lid for ultrasonic bath

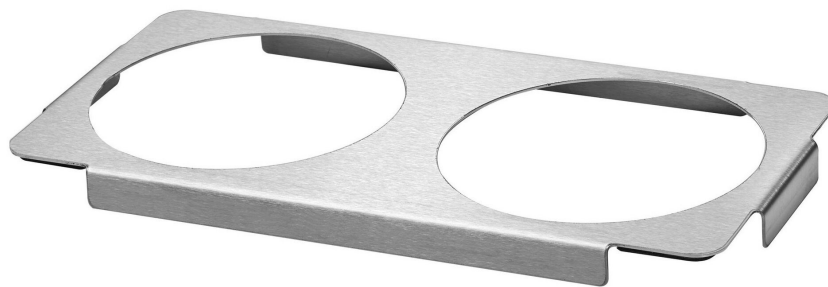

For the application of chemically aggressive liquids (e.g. TICKOPUR J 80 U), the removal of chemically aggressive or abrasive contaminations or the sonication of sample liquids, insert cups are used (indirect sonication). Because nothing may be placed on the bottom of the ultrasonic bath, BANDELIN offers suitable perforated lids/insert strips to hold the insert beakers. Further advantage: Several insert cups can be sonicated with different liquids at the same time (e.g. removal of dental cements with STAMMOPUR Z in one and cleaning with simultaneous chemical disinfection of dental instruments with STAMMOPUR DR 8 in a second insert cup).

- Stainless steel
- for the insertion of 2 insertion cups
- Suitable for RK 31/H

Item no.

WL10477

Model	DE 08
Manufacturer	BANDELIN
Manufacturer article no.	278
GPSR manufacturer data	BANDELIN electronic GmbH & Co. KG Heinrichstrasse 3-4 DE-12207 Berlin www.bandelin.com
Length	210 mm
Width	105 mm
Height	11 mm
Sales unit	1 piece
Content unit	1 piece
Incl. battery	no

