


## SONOREX hanging basket for ultrasonic bath, stainless steel, 200 x 110 x 40 mm

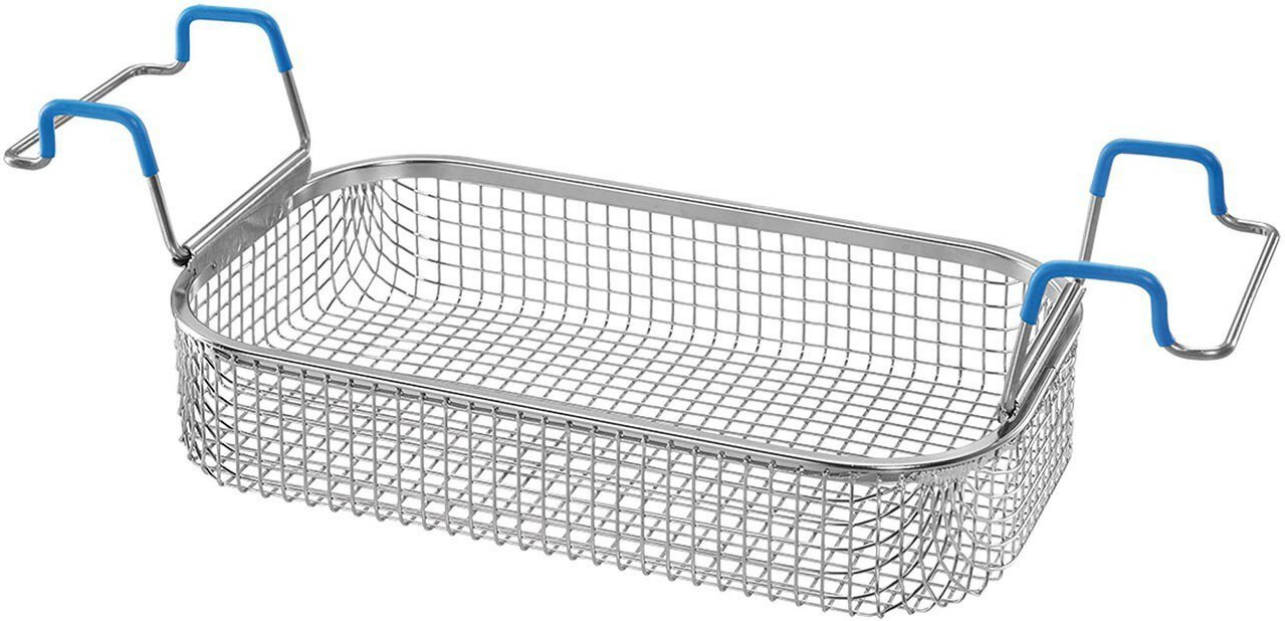

During ultrasonic cleaning, the items to be cleaned must not be placed on the bottom of the tank, as the vibrations of the oscillating systems may cause abrasion between the items to be cleaned and the tank itself. The use of a suitable cleaning basket is therefore recommended. The BANDELIN cleaning baskets are designed in such a way that the material to be cleaned is always in the optimum area of the cavitation field (1-2 cm from the floor) in the tank. With the aid of a basket, the items to be cleaned can also be conveniently placed in the unit and removed for rinsing after cleaning without the user coming into contact with the cleaning or disinfecting liquid.

- suitable for:
- DL: 102
- HDK: 102
- PDT: 100 /H, 102 H /H-RC
- RK: 100 /H /SH, 102 H /CH
- clear mesh width: 5 x 5 mm
- to hang in

|              |          |
|--------------|----------|
| Item no.     | WL10482  |
| Model        | K 3C     |
| Manufacturer | BANDELIN |

Products for the electronic industry



|                          |   |
|--------------------------|---|
| Manufacturer article no. | 3025  |
| GPSR manufacturer data   | BANDELIN electronic GmbH & Co. KG<br>Heinrichstrasse 3-4<br>DE-12207 Berlin<br>www.bandelin.com |
| Length                   | 200 mm  |
| Width                    | 110 mm  |
| Height                   | 40 mm   |
| Sales unit               | 1 piece   |
| Content unit             | 1 piece   |

