


SD 05 Insert cup cpl. For ultrasonic bath

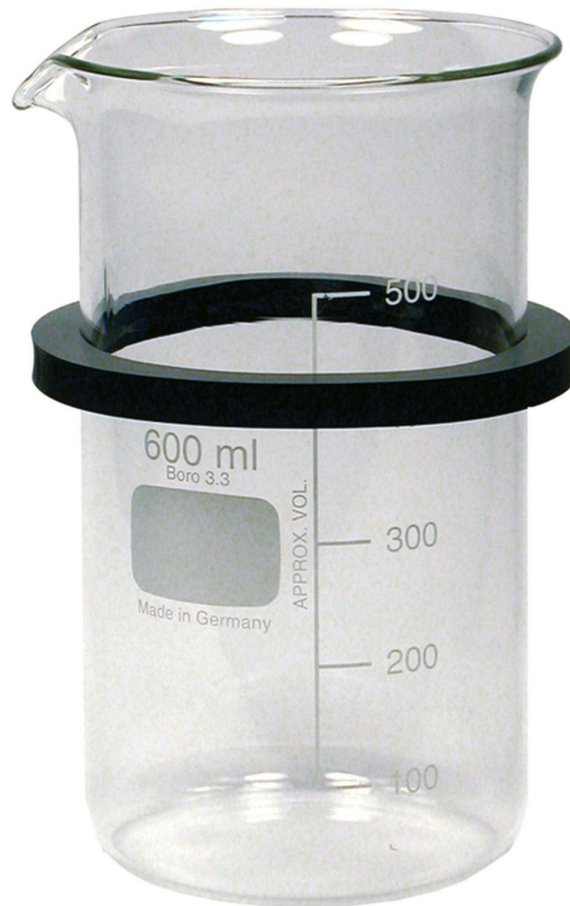

Glass, 600 ml, for insertion into perforated lid/insertion strip with ring GR 04. Insert beakers are used for the application of chemically aggressive liquids (e.g. TICKOPUR J 80 U), the removal of chemically aggressive or abrasive contaminants or the sonication of sample liquids (indirect sonication). Because nothing may be placed on the bottom of the ultrasonic bath, BANDELIN offers suitable perforated lids/insert strips to hold the insert beakers. Further advantage: Several insert cups can be sonicated with different liquids at the same time (e.g. removal of dental cements with STAMMOPUR Z in one and cleaning with simultaneous chemical disinfection of dental instruments with STAMMOPUR DR 8 in a second insert cup). Please note that in order to transfer the ultrasound to the insert cup, the surface tension of the water must be reduced and therefore a cleaning solution (e.g. TICKOPUR R 33 3%) must be used as coupling liquid / contact liquid.

□ for suspension in perforated cover DE08

□ for RK 31 / RK 31 H

□ D = 76 mm, H = 110 mm

□ Capacity = 400 ml

- Glass, complete with retaining ring
- for indirect cleaning of small parts

Item no.	WL10498
Model	SD 05
Manufacturer	BANDELIN
Manufacturer article no.	575
GPSR manufacturer data	BANDELIN electronic GmbH & Co. KG Heinrichstrasse 3-4 DE-12207 Berlin www.bandelin.com
Length	76 mm
Width	80 mm
Height	150 mm
Sales unit	1 piece
Content unit	1 piece

