


SONOREX hanging basket for ultrasonic bath, stainless steel, 545 x 450 x 50 mm

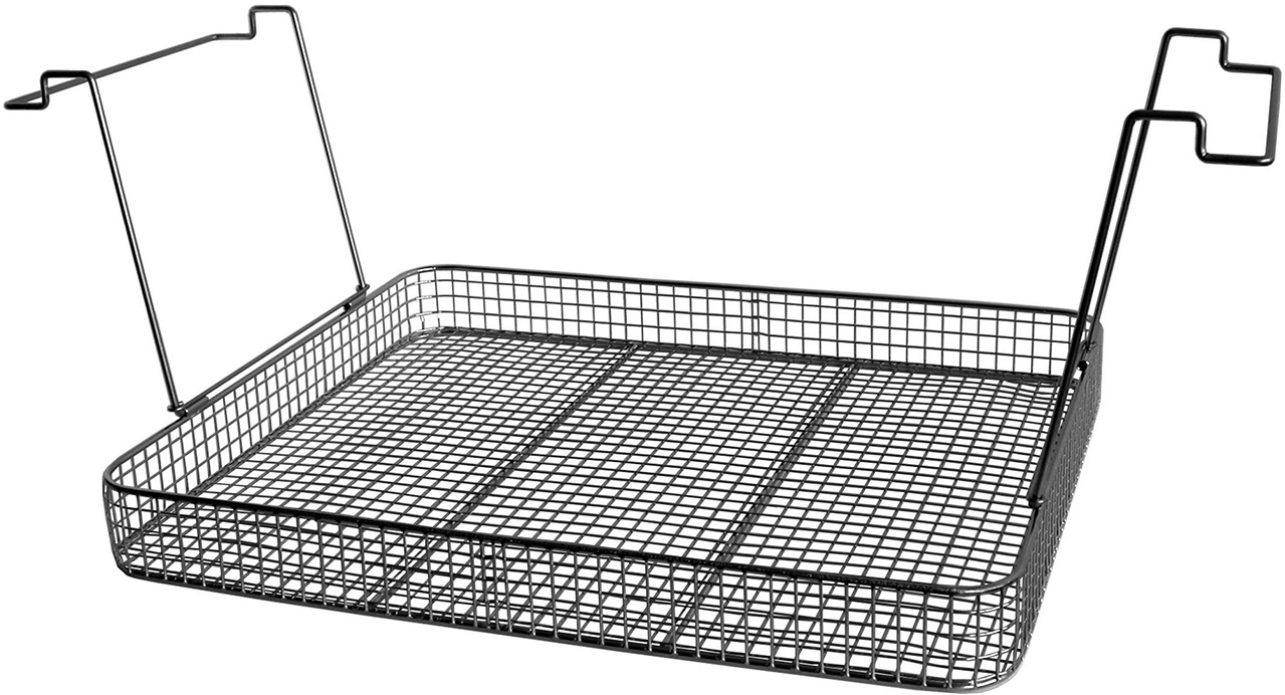

During ultrasonic cleaning, the items to be cleaned must not be placed on the bottom of the tank, as the vibrations of the oscillating systems may cause abrasion between the items to be cleaned and the tank itself. The use of a suitable cleaning basket is therefore recommended. The BANDELIN cleaning baskets are designed in such a way that the material to be cleaned is always in the optimum area of the cavitation field (1-2 cm from the floor) in the tank. With the aid of a basket, the items to be cleaned can also be conveniently placed in the unit and removed for rinsing after cleaning without the user coming into contact with the cleaning or disinfecting liquid.

- suitable for:
- DT: 1050 CH
- RK: 1050 CH
- Mesh size: 5 x 5 mm

item number	WL19127
model	K 50C
manufacturer	BANDELIN
manufacturer item number	138

Products for the electronic industry



length	545 mm
width	450 mm
height	50 mm
order unit	1 piece
content unit	1 piece