


SONOREX hanging basket for ultrasonic bath, stainless steel, $460 \times 100 \times 50 \text{ mm}$

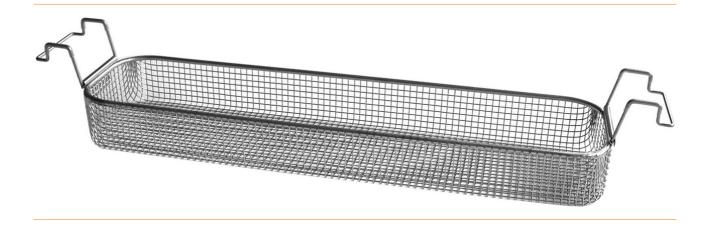

During ultrasonic cleaning, the items to be cleaned must not be placed on the bottom of the tank, as the vibrations of the oscillating systems may cause abrasion between the items to be cleaned and the tank itself. The use of a suitable cleaning basket is therefore recommended. The BANDELIN cleaning baskets are designed in such a way that the material to be cleaned is always in the optimum area of the cavitation field (1-2 cm from the floor) in the tank. With the aid of a basket, the items to be cleaned can also be conveniently placed in the unit and removed for rinsing after cleaning without the user coming into contact with the cleaning or disinfecting liquid.

suitable for: DT: 156, RK: 156clear mesh width: 5 x 5 mm

• to hang in

item number	WL27058
model	K 6L
manufacturer	BANDELIN
manufacturer item number	202
length	460 mm
width	100 mm
height	50 mm
order unit	1 piece
content unit	1 piece









24.04.2024